



Community Plus™

AHSM-072 Rev 2.0 (25.10.24)


Static Replacement Mattress



The **Community Plus** mattress has been designed for users considered to be at '**High Risk**' of developing pressure ulcers and developed to meet the demands of the community and care home environments.

Comprising of a high specification castellated reflex foam, the **Community Plus** mattress is cost-effective and durable, delivering exceptional levels of pressure reduction, comfort and support.

- Castellated reflex foam offers excellent pressure reduction, comfort, and support
- Internal air channels in the foam's surface help to reduce the buildup of heat and moisture on the skin
- Smooth multi-stretch PU cover is both waterproof and vapour permeable, reducing shear and friction while reducing the risk of fluid ingress, supporting infection control protocols
- The non-turn mattress system reduces manual handling risks
- Suitable for electric profiling bed frames
- Available in single and double sizes, bespoke sizes are also available upon request

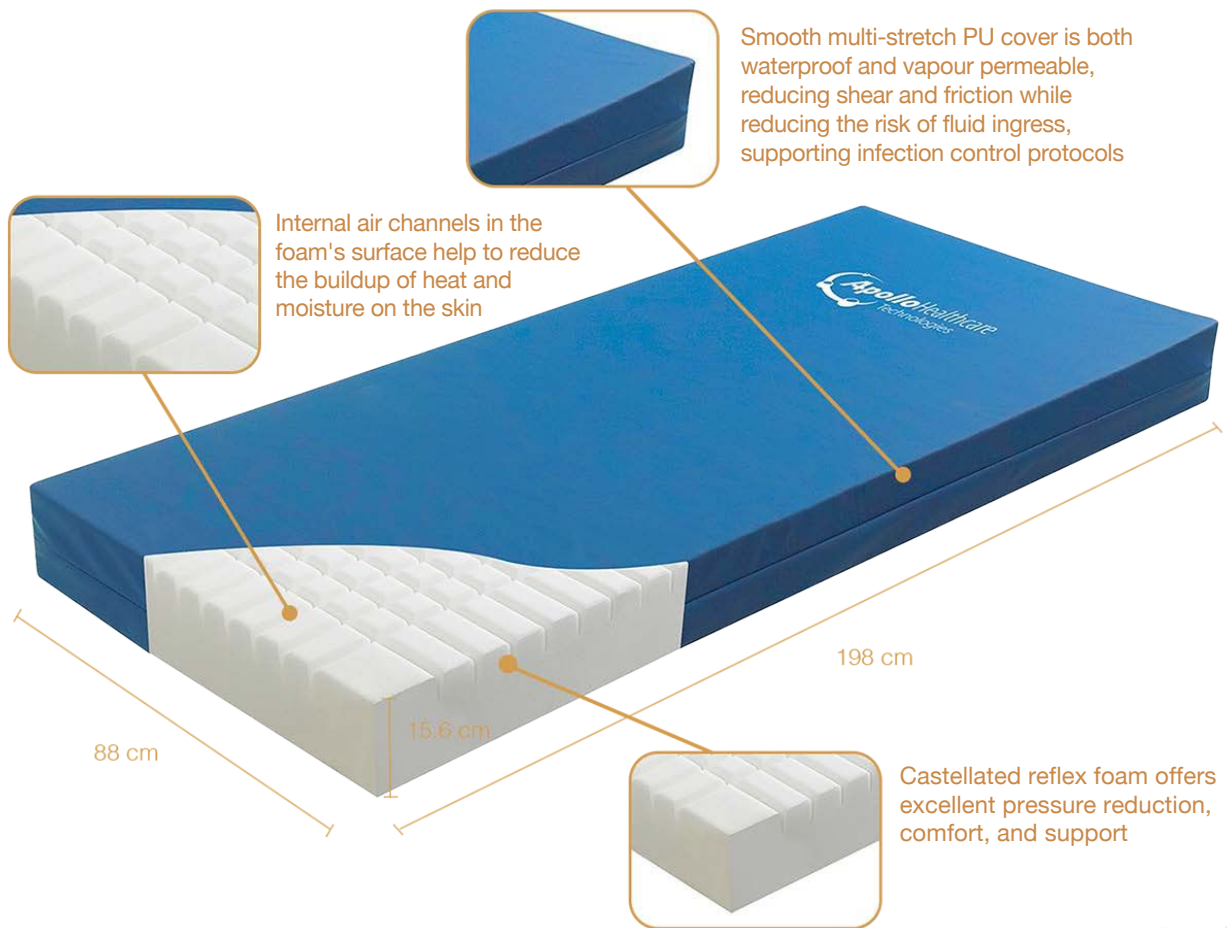

Specifications

Product Code - APH003

Dimensions: 198 x 88 x 15.6 cm

Mattress Weight: 12 kg

Maximum Weight Limit: 36 stone / 229 kg

Risk Factor: High Risk

Warranty: 2 years

Safety Standards: BS EN 597-1 Cig, BS EN 597-2 Match, BS7177 Ignition Source 5

Due to ongoing research and development, Apollo Healthcare Technologies Ltd, reserve the right to change specifications without prior notice.

If you require more information, please call our sales team on Tel: +44 (0) 1924 614567



Apollo Healthcare Technologies Limited, Unit A1, Gildersome Spur, Leeds, West Yorkshire, LS27 7JZ
Tel: +44 (0) 1924 614567 Fax: +44 (0) 1924 607480 Email: sales@apollo-ht.co.uk www.apollo-ht.co.uk